
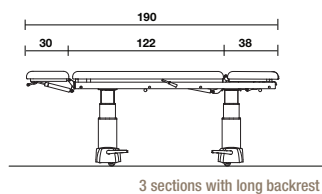
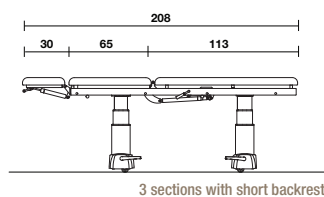
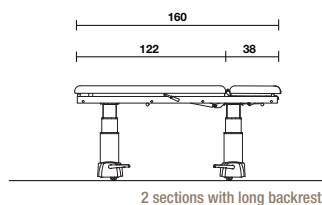


Helio Examination table

15064 The examination table consists of a mattress platform with two or three sections. The open undercarriage increases the ergonomic ease for both the patient and the doctor.

Open undercarriage

- High/low adjustment by means of 2 telescopic lifting columns.
- Stable, absolute vertical high/low adjustment.
- The castors are hidden by means of a synthetic cover cap with a grey uni colour.
- Twin castors Ø 50 mm, braked 2 per 2.
- Design pedal to activate the brake mechanism.
- The design brake pedal can be activated on both sides of the bed.
- The open carriage is very accessible and easy to clean.

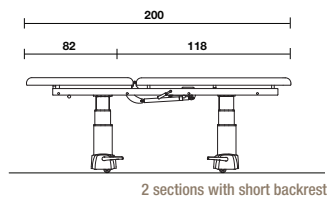

Mattress platform

- Mattress platform: choice between 2 or 3 sections.
- The mattress platform is upholstered with fire-retardant foam. The foam is 40 mm thick and has a density of 45 kg/m³.
- Finished with imitation leather.
- Two-section mattress platform with short backrest or long backrest.
- Three-section mattress platform with headrest and choice between short or long backrest.
- The backrest can be manually adjusted up to ± 85° by means of a gas spring.
- The headrest can be inclined manually upwards or downwards by means of a gas spring (-40° to +40°). (only with 3-sections)
- The seat support can be manually adjusted by means of a telescopic adjust.

Gas spring

- The headrest and backrest can be inclined by means of a gas spring.
- The gas spring that inclines the headrest can be activated by means of a lever that is situated on the side of the headrest.
- The backrest is inclined by means of a lever situated on both sides of the examination table.

Electric manipulations

- The high/low adjustment can be activated by means of a foot pedal that can be positioned freely on the floor.
- The electric backrest adjustment (*) can be activated by means of a double foot pedal.

Colours

We would like to refer to our extensive colour gamut. Sample card + samples available on request.

Accessories

We would like to refer to our accessory list.



Dimensions and prestations

- Width: 70 cm
- Length: depends on the type of patient platform
- Height variation: ± 52 cm - ± 92 cm
- Weight (without options): 74 kg
- Safe working load: 250 kg
- Headrest (*) can be inclined between -40° and +40°
- Backrest (*) can be adjusted up to 85°

Materials and finishing

- Metal: epoxy-coating RAL 7035.
- Fire-retardant foam.
- Upholstery: see collection.
- Withstands the common cleaning products and disinfectants.

For further information about the **used materials, constructions and maintenance**: please consult our material files.

Options

- Gynaecology model (only with long backrest).
- Handles (only with long backrest).
- DIN rail (foot, head, head and foot).
- IV rod (fixation included).
- Armrest (fixation included).
- Medical lamp.
- Sideguards.
- Paper roll holder.
- Electric adjustment of the backrest.
- Opening for the nose in the headrest or in the 2-sections table with short backrest.
- Neck cushion and/or lumbar support cushion.
- Equipotential bonding.
- Comfort Plus mattress platform.
- Battery (Li-Ion ou lead acid).
- Removable tray (only with long backrest).
- Transparent cover for footrest (not available on versions with long backrest or Comfort Plus mattress platform).
- Trendelenburg position by means of a handset.

(*) Option

Order here!

PRAXISDIENST
Medical Supplies since 1953

